

MKS® Body Belt plus

☎ 00800 03032900

Lumbar Bandage

Sizes



This product is available in pediatric sizes.

Indications

- Lower back pain
- Lumbago
- Spondylarthrosis
- Sacroiliac joint irritation
- Lumbar muscle insufficiency

Features

- Semi-elastic bandaging material
- Lateral reinforcing stays
- Firm abdominal piece with hook and loop fastener
- Back pocket suitable for silicone pad stimulating blood circulation or MKS® Osteo back plate
- Additional tensioning straps



Mode of Operation


Due to its semi-elastic bandaging material the MKS® Body Belt plus adjusts comfortably to the body. Through building up circular pressure the lumbar vertebral column is supported.

The firm fastener element provides the necessary abdominal pressure which can be enhanced through the additional tensioning straps. Optionally, blood circulation can be stimulated by inserting the silicone pad into the integrated back pocket.

Size chart and article numbers

Size	Body circumference	Abdominal height	Height of back piece	Art.-No.	Reduced height		
					Abdominal height	Height of back piece	Art.-No.
1	60 - 70 cm	18 cm	22 cm	582-201-S	14 cm	18 cm	582-201N-S
2	70 - 80 cm	18 cm	23 cm	582-202-S	14 cm	19 cm	582-202N-S
3	80 - 90 cm	18 cm	23 cm	582-203-S	14 cm	19 cm	582-203N-S
4	90 - 105 cm	18 cm	24 cm	582-204-S	14 cm	20 cm	582-204N-S
5	105 - 120 cm	18 cm	24 cm	582-205-S	14 cm	20 cm	582-205N-S
6	120 - 135 cm	18 cm	24 cm	582-206-S	14 cm	20 cm	582-206N-S

Silicone pad optionally available. (Art.-No.: 350-014-S)

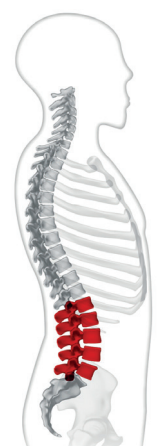

L1 - L5



Spinal
orthoses

Other sizes available on request.

www.albrechtgmbh.com

