

MKS® Lumbago plus

Lumbar Bandage with silicone pad stimulating blood circulation

Indications

- Severe lower back pain
- Severe lumbago
- Relapsing lumbago
- Osteochondrosis
- Mild lumbar sciatica
- Mild radicular syndrome
- Spondylolysis without displacement of vertebrae
- Mild lumbar deformity (with displacement of vertebrae)
- Spondylarthrosis
- Sacroiliac joint irritation
- Lumbar muscle insufficiency

Features

- Semi-elastic bandaging material
- Lateral reinforcing stays
- Firm abdominal piece with hook and loop fastener
- Back pocket suitable for silicone pad stimulating blood circulation or MKS® Osteo back plate
- Additional tensioning straps
- Silicone pad stimulating blood circulation




Mode of Operation

Due to its semi-elastic bandaging material the MKS® Lumbago plus adjusts comfortably to the body. Through building up circular pressure the lumbar vertebral column is supported. The firm fastener element provides the necessary abdominal pressure which can be enhanced through the additional tensioning straps. The inserted silicone pad in the integrated back pocket stimulates blood circulation.



Silicone Pad

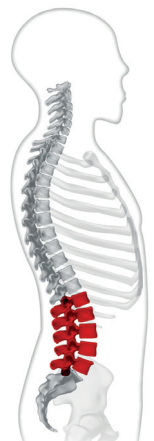

L1 - L5



Size chart and article numbers

Size	Body circumference	Abdominal height	Height of back piece	Art.-No.
1	60 - 70 cm	18 cm	22 cm	582-501-S
2	70 - 80 cm	18 cm	23 cm	582-502-S
3	80 - 90 cm	18 cm	23 cm	582-503-S
4	90 - 105 cm	18 cm	24 cm	582-504-S
5	105 - 120 cm	18 cm	24 cm	582-505-S
6	120 - 135 cm	18 cm	24 cm	582-506-S

