

Lumbar bandage with firm bridging frame

Indications

- Severe lumbar sciatica with muscular failures in severe intervertebral disc protrusion/prolapse without absolute indication for surgery
- Severe radicular lumbar syndrome, pseudo-radicular lumbar syndrome
- Spondylolithesis with relapsing sciaticae
- Severe lumbar deformity caused by facet syndrome/arthrosis
- Condition after intervertebral disc surgery without residual paresis in cleared intervertebral disc, decompression of the spinal canal one-layer without deformity, intervertebral disc surgery one-layer
- Spinal stenosis without strong paresis after short walking distance

Features

- Semi-elastic bandaging material
- Firm abdominal piece with hook and loop fastener
- Bridging rods (spring band steel)
- Lateral reinforcing stays
- Additional lateral double straps
- Firm bridging frame




Mode of Operation

The MKS® Pontsana® stabil functions according to the bridging principle. The wide bandage with the additional lateral double straps and the firm bridging frame straightens and hence relieves the lumbar vertebral column.

Size chart and article numbers

Size	Body-circumference	Abdominal height	Height of back piece	Art.-No.	Reduced height		
					Abdominal height	Height of back piece	Art.-No.
1	60 - 70 cm	18 cm	32 cm	582-881-S	16 cm	26 cm	582-881N-S
2	70 - 80 cm	18 cm	32 cm	582-882-S	16 cm	26 cm	582-882N-S
3	80 - 90 cm	18 cm	34 cm	582-883-S	16 cm	26 cm	582-883N-S
4	90 - 105 cm	18 cm	34 cm	582-884-S	16 cm	26 cm	582-884N-S
5	105 - 120 cm	18 cm	38 cm	582-885-S	16 cm	26 cm	582-885N-S
6	120 - 135 cm	18 cm	38 cm	582-886-S	16 cm	26 cm	582-886N-S
51	105 - 120 cm	18 cm	34 cm (reduced)	582-8851-S			
61	120 - 135 cm	18 cm	34 cm (reduced)	582-8861-S			

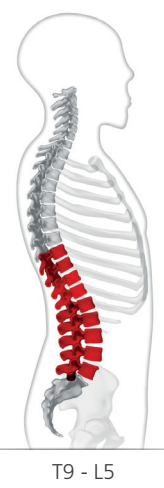

Spinal
orthoses

Other sizes available on request.

